


XP SUPPORT OVERLAY

Overlays are an effective solution for a wide variety of conditions around client comfort, pressure care and sleep climate.

FIRM FEEL


LOW RISK

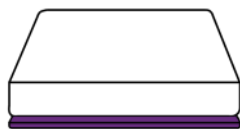


Use as an underlayer. This prevents bottom out and adds support.



Use to make an existing mattress firmer. This can help with mobility and transfers.

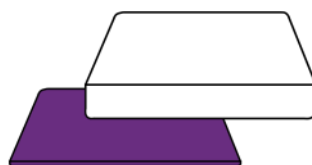

Use as a crash mat. Can be placed beside bed to reduce risk of injury.

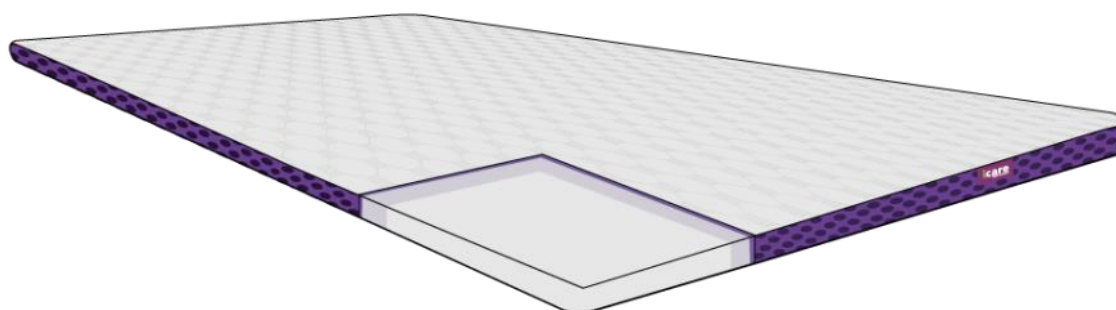


Solid Support Foam Base

XP Support foam is a high-density open cell material that breathes to allow air circulation but is high density to prevent bottoming out. This material will easily contour to any profiling bed.



Ventilated Base Fabric

Vant fabric is a three-dimensional polypropylene material that has high air permeability and is anti-bacterial. We use this on the side wall of this overlay cover to greatly increase air movement.

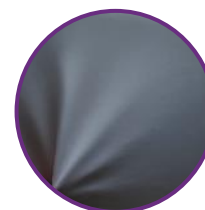

Medical Covers Available

Zerotec fabric is a medical barrier fabric that's has high elasticity in every direction which reduces shear on pressure redistribution surfaces. It is waterproof and wipeable but has microporous breathable properties.

Alternating Air Mattresses

Overlays are frequently used with an alternating air mattress. When an air mattress is used on a profiling bed, the air pressure can change when the bed positions are adjusted. The XP Support Overlay is suitable to use under an alternating air mattress. When an air mattress loses pressure in profiled positions, patients are not left lying on the bars or mesh of beds. This can prevent pressure sores or injuries.

CODES & SIZES:



SINGLE

Code: ICXOS
Mattress Size: 1900 x 900



LONG SINGLE

Code: ICXOLS
Mattress Size: 2030 x 900



KING SINGLE

Code: ICXOKS
Mattress Size: 2030 x 1070



DOUBLE

Code: ICXOD
Mattress Size: 1900 x 1350



QUEEN

Code: ICXOQ
Mattress Size: 2030 x 1520

KEY INFORMATION:



SWL Minimum 50kg
SWL Maximum 250kg



5 Year Warranty



50mm Depth


*Measurements in Millimetres

ACTIVEX™ OVERLAY

The Icare ActiveX™ Pressure Care Overlay is designed to provide added pressure relief and comfort on an existing mattress. This overlay is heat and pressure sensitive causing them to soften around pressure points and inflamed areas.



Use to increase comfort on mattresses especially for short-term illness, recovery from surgery or for added comfort.

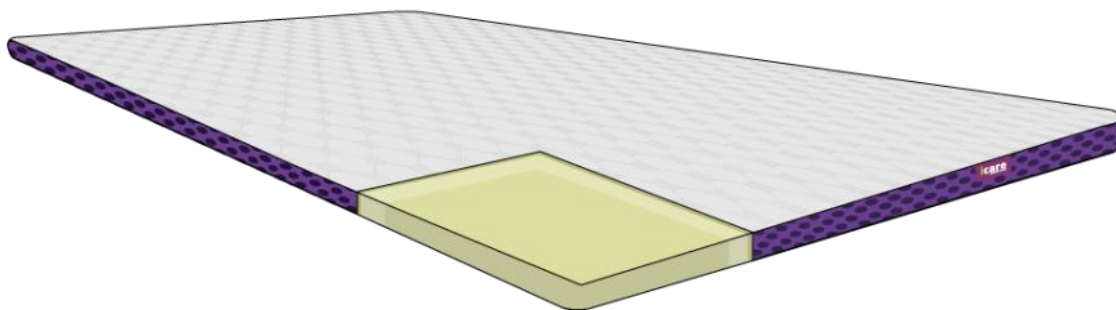

Alternating Air Mattresses

Safe to use on alternating air mattresses. ActiveX™ Overlay's are often used to increase comfort and reduce heat loss when used with alternating air mattresses.



Pressure Redistribution

The overlay features the ActiveX™ 84 material which is a high grade of elastic type foam. ActiveX™ is responsive to temperature around inflamed areas and will become softer as it becomes warmer. So around a pressure point it will always soften and allow the blood to circulate freely.



Ventilated Base Fabric

Vant fabric is a three-dimensional polypropylene material that has high air permeability and is anti-bacterial. We use this on the side wall of this overlay cover to greatly increase air movement.



Medical Covers Available

Zerotec fabric is a medical barrier fabric that's has high elasticity in every direction which reduces shear on pressure redistribution surfaces. It is waterproof and wipeable but has microporous breathable properties.

CODES & SIZES:



SINGLE

Code: ICTS
Mattress Size: 1900 x 900



LONG SINGLE

Code: ICTLS
Mattress Size: 2030 x 900



KING SINGLE

Code: ICTKS
Mattress Size: 2030 x 1070



DOUBLE

Code: ICTD
Mattress Size: 1900 x 1350



QUEEN

Code: ICTQ
Mattress Size: 2030 x 1520

*Measurements in Millimetres

SOFT FEEL

HIGH RISK

This rating is based on our average peak pressure of <22mmHg. For some high risk clients Icare Zerotec covers are recommended.

Icare risk category ratings are based on average peak pressure, Therapist recommendation, case studies, industry comparison and the Waterlow scale. For some high risk clients Icare Zerotec covers are recommended. The best way to assess a client on a pressure reducing surface is to monitor pressure points with visual checks over a period of 5 days.

KEY INFORMATION:



SWL Minimum 30kg
SWL Maximum 215kg



5 Year Warranty



50mm Depth

ACTIVEFLOAT™ OVERLAY

Our Icare ActiveFloat™ Overlay is the latest technology is an extremely elastic and highly breathable material designed to reduce shear and offer pressure relief.

SOFT FEEL


HIGH RISK

This rating is based on our average peak pressure of <22mmHg.


Icare risk category ratings are based on average peak pressure, Therapist recommendation, case studies, industry comparison and the Waterlow scale. The best way to assess a client on a pressure reducing surface is to monitor pressure points with visual checks over a period of 5 days.



Premium Overlay that reduces shear in every direction with increased user comfort.

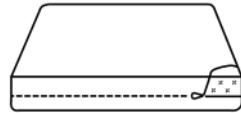

SlideSkin™ reduces shear enables the cover to move or slide with the user movements.

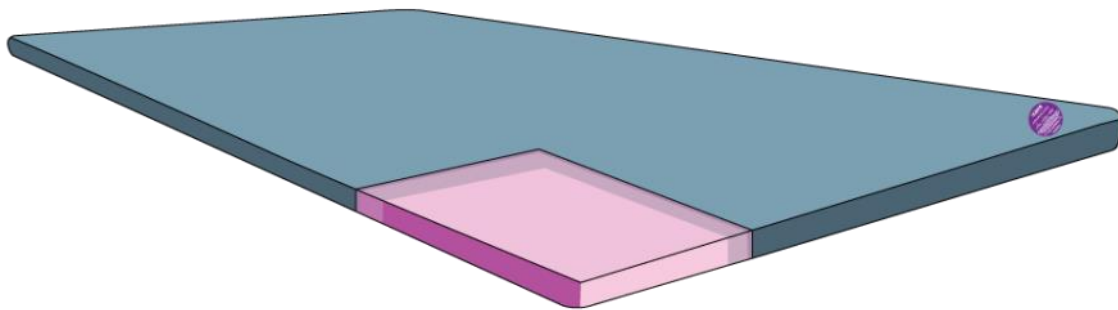

Pressure Redistribution

The overlay combines the features the new ActiveFloat™ foam which reduces shear whilst remaining fully pliant under maximum compression. This provides a product that is responsive to temperature around inflamed areas and will become softer as it becomes warmer. So around a pressure point it will always soften and allow the blood to circulate freely.



SlideSkin™ Anti-Shear

The unique composition of SlideSkin™ technology reduces shear with a second membrane. This enables the cover to move or slide with the users movements.



Medical Cover

Zerotec fabric is a medical barrier fabric that's has high elasticity in every direction which reduces shear on pressure redistribution surfaces. It is waterproof and wipeable but has microporous breathable properties. This cover has welded seams and waterfall flaps for fluid control.

CODES & SIZES:



SINGLE

Code: ICOFS
Mattress Size: 1900 x 900



LONG SINGLE

Code: ICOFLS
Mattress Size: 2030 x 900



KING SINGLE

Code: ICOFKS
Mattress Size: 2030 x 1070



DOUBLE

Code: ICOFD
Mattress Size: 1900 x 1350



QUEEN

Code: ICOFQ
Mattress Size: 2030 x 1520

*Measurements in Millimetres

KEY INFORMATION:



SWL Minimum 30kg
SWL Maximum 215kg



5 Year Warranty



50mm Depth